TD.80396 Iron Deficient Diet

Formula	g/Kg
Casein, low Cu & Fe	200.0
DL-Methionine	3.0
Sucrose	549.99
Corn Starch	150.0
Corn Oil	50.0
Mineral Mix, Fe Deficient (81062)	35.0
Vitamin Mix, AIN-76A (40077)	10.0
Choline Bitartrate	2.0
Ethoxyquin, antioxidant	0.01

Key Features

- Purified Diet
- AIN-76A Modification
- Iron Deficient

Key Planning Information

- Products are made fresh to order
- Store product at 4°C or lower
- Use within 6 months (applicable to most diets)
- Box labeled with product name, manufacturing date, and lot number
- Replace diet at minimum once per week
 More frequent replacement may be advised
- Lead time:
 - · 2 weeks non-irradiated
 - · 4 weeks irradiated



Product Specific Information

- 1/2" Pellet or Powder (free flowing)
- Minimum order 3 Kg
- Irradiation not advised
 - · Contact a nutritionist for recommendations

Footnote

The diet contains approximately 2-6 ppm Fe. To limit background iron, cellulose is omitted, and reagent grade, pretested calcium phosphate is used in the mineral mix. TD.80394 is a possible ingredient matched control with ~48 ppm added iron from ferric citrate. For questions about this or other iron adjusted diets askanutritionist@harlan.com.

Selected Nutrient Information¹

	% by weight	% kcal from
Protein	17.7	17.8
Carbohydrate	69.9	70.4
Fat	5.2	11.8


Kcal/a 4.0

¹Values are calculated from ingredient analysis or manufacturer data

Teklad Diets are designed & manufactured for research purposes only.

Speak With A Nutritionist

- (800) 483-5523 option 4
- askanutritionist@harlan.com

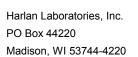





www.harlan.com

Options (Fees Will Apply)

- Rush order (pending availability)
- Irradiation (see Product Specific Information)
- Vacuum packaging (1 and 2 Kg)

Contact Us

- · Obtain Pricing · Check Order Status ·
- tekladinfo@harlan.com
- (800) 483-5523 option 3

International Inquiry (Outside USA or Canada)

• askanutritionist@harlan.com

Place Your Order (USA & Canada)

- · Please Choose One ·
- https://webapps.harlan.com/teklad-orders/
- tekladinfo@harlan.com
- (800) 483-5523 option 3
- (608) 277-2066 facsimile

Helping you do research better